cosmobelle

cavitation

your silhouette is up to you

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your silhouette is up to you

USER MANUAL

Congratulations, First of all we would like to thank you for purchasing SILUET W&M and for the trust you have placed in Cosmobelle. By means of this handy User Manual and following its simple instructions you will achieve results comparable to those obtained from the hands of a professional beautician.

We strongly recommend that you read this User Manual thoroughly, reading and understanding it will help you use SILUET W&M effectively and safely.

CONTENTS

CHAPTER 1

- contents of the box
- SILUET W&M assembly

CHAPTER 2

- basic features of SILUET W&M
- operating instructions

CHAPTER 3

- application modes

CHAPTER 4

- correct, safe use of SILUET W&M

CHAPTER 5

- frequently asked questions about SILUET W&M

CHAPTER 6


- important safety instructions
- parts list
- certificates
- warranty
- technical specifications

CHAPTER 1

Contents of the box

Inside this box you will find everything you need to enjoy SILUET W&M from now onward. Carefully check that the pack contains all of the elements on the following list.

- User manual
- DVD
- Main unit
- Hand unit
- Mains cable
- Cavitation gel
- Post cavitation Lotion
- Toning serum
- Anti-cellulite serum

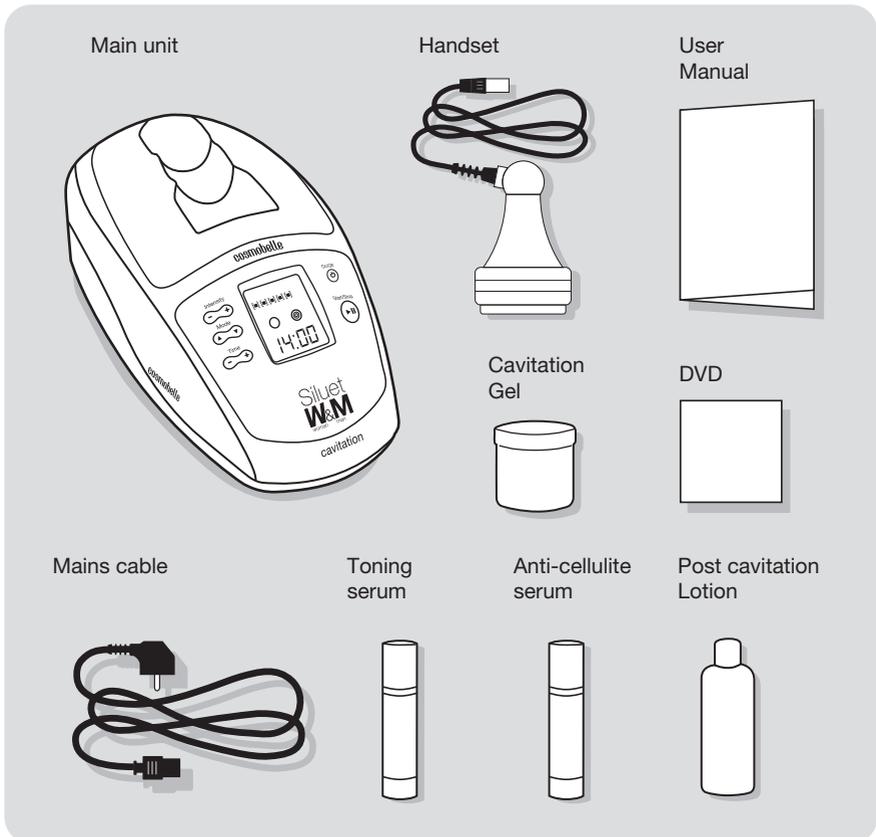CHAPTER 2

Basic features of Siluet W&M

SILUET W&M is a hi-tech device that allows us to selectively and efficiently remove the build-up of fat (cellulite, orange skin syndrome) from our body.

The technology used is that of ultrasound emission which produces the physical phenomenon known as cavitation, causing the removal of fat cells.

This technology with proven results is widely used in beauty salons all over the world in order to fight cellulite painlessly and non-invasively.

En

Characteristics of the SILUET W&M model

- Time: Programmable from 1 to 15 minutes.
- Intensity: Programmable at 6 intensity levels.
- Mode: 1 pulsed mode and 1 continuous mode.

With this unit you only have to worry about programming in three simple functions: intensity, time and mode.

The intensity function controls how deep the sound waves penetrate your skin. The greater the intensity, the deeper the waves reach.

The time function allows you to choose how long you wish to work with SILUET W&M. For reasons of safety the time function allows you to choose up to 15 minutes, which corresponds to the maximum recommended time for working a specific area. However, you can reduce this maximum time by 1 minute intervals. Once the selected time has elapsed, SILUET W&M will inform you with a beep.




SILUET W&M also lets you choose the work mode, helping you to select the time interval during which the waves reach your skin. In pulsed mode, the waves are emitted intermittently, whereas in continuous mode the emission is continuous, which means it works harder.

Siluet W&M operating instructions

Precautions: Consult your doctor before beginning to use this appliance. This warning is especially important for people with health problems. Use by those under 18 years of age is not recommended. Keep this appliance out of children's reach.

Starting

- Connect the hand unit to the main unit. To do this, take the plug for the handset and insert into the socket on the main unit as shown in fig. 2.1
- Plug one end of the mains cable into the socket on the machine and the other end into the electric wall socket (220V). fig. 2.2
- Press the switch at the back of the unit to start it (position I). A red light below the screen will light up. fig. 2.3
- Press the  On/Off button to switch the machine on. The screen will switch on and green light below it.

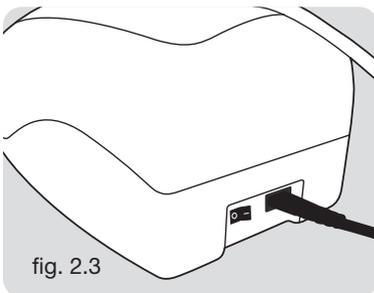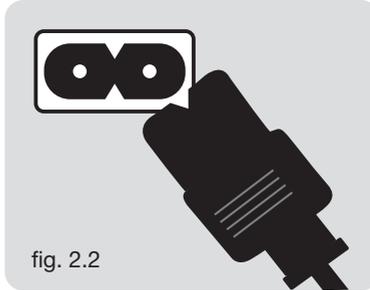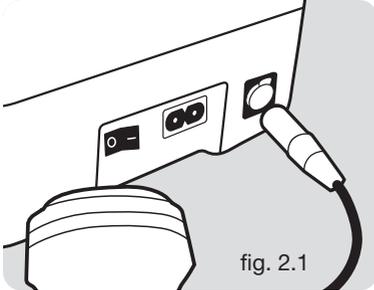Monitor elements

The SILUET W&M monitor consists of a display screen and five buttons or keys.

- The display provides information about the cavitation session.
- The keys have the following functions:
 - **+** and **- Intensity** (change intensity)
 - **+** and **- Time** (Time settings)
 - **▲** and **▼ Mode** (change mode)
 - **⏻ On/Off**
 - **▶ || Start/Stop**

Cavitation session

There are 3 settings for the massage:

- **Time:** measures the duration of the session in minutes and seconds. You can select between 1 to 15 minutes (in steps of 1 minute).
- **Intensity:** establishes the intensity level for the session. SILUET W&M lets you choose from among 6 levels.
- **Mode:** Pulse type (impulse emission frequency). SILUET W&M allows you to choose between pulsed and continuous mode.

Beginning a session

Once you have plugged the machine into the mains supply and switched it to **position I** press the **On/Off** button and the following default settings will appear on the monitor:

- Time: 10 minutes (10:00:00)
- Intensity: Level 4 **[■][■][■][■]**
- Mode: Continuous **○**

From here you can modify the time, intensity and mode settings to establish the type of session you desire by doing the following:

Changing session duration

Once you are on the default setting screen you can modify the duration of the session by using the **+** and **– (Time)** keys.

The displayed default value is 10:00:00 (10 minutes) but this can be changed in steps of 1 minute between the values 1:00 and 15:00.

Changing massage intensity

Once you have selected the desired time, you can then set the intensity. To do this press the **+** and **– (Intensity)** keys.

The intensity level that appears by default, or from the last session, can be changed by selecting between levels 1 and 6 (1 to 6 squares). The intensity setting can be changed during operation by pressing the **+** and **–** keys.

Changing work mode

You can use the **▲** and **▼** keys (**Mode**) to change the work mode for your session.

By default, the machine operates in **continuous** mode, which is identified with the flashing empty circle located on the left side of the screen . **Pulsed** mode, is identified with the concentric circle icon located on the right side of the screen .

Finally, grip the handset firmly and place it onto the area that you wish to work, then press the **On/Off**  button to start the session. Massage as shown on the DVD or in the “application modes” section.

If you wish to stop the session, just press the **▶ || Start/Stop** button and this will stop the session, resetting the counter to the previously selected time.

When the time counter reaches **0:00** the unit will stop and emit a beep to tell you that the selected session has finished. The counter will return to the previously selected time.

CHAPTER 3

Siluet W&M application modes

Always apply an abundant amount of gel to the area that you wish to work and begin to treat it by sliding the handset slowly and without pressing, as shown in the diagrams.

If you need to purchase gel or any of our pre-cavitation creams you can find them at the retail outlet where you purchased SILUET W&M or find other sales outlets at the website: www.cosmobelle.com

In each of the following sections you will find a table with application protocols. Before beginning to use the appliance, you must determine the level of fat in each area to be treated.


En

- A** Large accumulation of high density fat
- B** Large accumulation of medium density fat
- C** Medium accumulation of medium density fat
- D** Small accumulation of fat
- E** Orange peel syndrome
- F** Maintenance session

Central abdomen

Remember, never apply SILUET W&M to the navel, breasts, the thorax area or the genitals.

Apply an abundant amount of gel on your abdomen and begin to massage the area by sliding the handset slowly and without pressing, as shown in the diagrams.

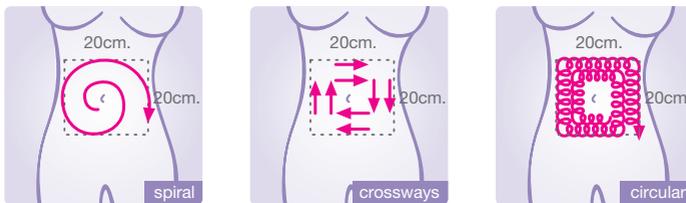
For best results we recommend applying SILUET W&M by combining different massage techniques (spiral, crossways, circular) during different sessions.

We recommend lying down for a more comfortable SILUET W&M application.

Level	Intensity	Time	Mode
A	4	10'	☉
B	4	10'	○
C	3	10'	○
D	2	8'	○
E	2	8'	○
F	1	8'	○

Sides

Apply gel over the area to be worked using SILUET W&M.

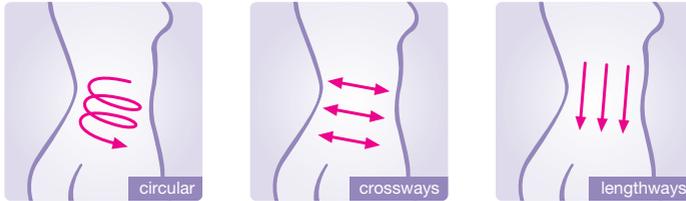
Circular and longitudinal massage movements should be applied from top to bottom as shown in the diagrams. It is important to massage in the correct direction as this will improve drainage. Combine the various massage techniques during different sessions.

Level	Intensity	Time	Mode
A	5	10'	○
B	5	10'	⊙
C	4	10'	⊗
D	2	8'	⊙
E	2	5'	⊙
F	1	5'	⊙

Hips

Apply gel over the area to be worked using SILUET W&M.

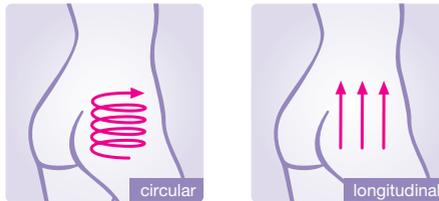
Circular and longitudinal massage movements should be applied from the bottom upward as shown in the diagrams. During treatment concentrate on the area with the most fat.

Level	Intensity	Time	Mode
A	4	8'	⊙
B	4	8'	⊙
C	3	8'	⊙
D	3	8'	⊙
E	2	5'	⊙
F	1	5'	⊙

Buttocks

Apply gel over the area to be worked using SILUET W&M.

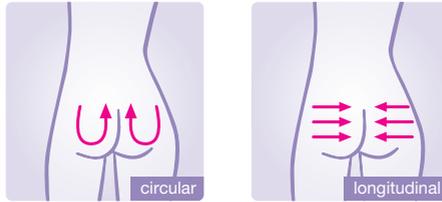En


Longitudinal and circular massage should be done in the direction of the arrows shown in the diagrams.

Level	Intensity	Time	Mode
A	6	10'	○
B	6	10'	○
C	5	8'	⊙
D	4	8'	⊙
E	3	10'	⊙
F	1	10'	⊙

Inner thigh

Remember, never apply SILUET W&M from the knee to the foot, there is no significant accumulation of fat in this area.

Apply gel over the area to be worked using SILUET W&M.

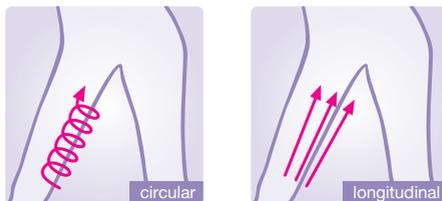
Both the longitudinal and circular massage should be done from the bottom upward as shown in the diagrams. Follow these simple instructions to help you get the best results.

Level	Intensity	Time	Mode
A	4	10'	⊙
B	4	10'	⊙
C	3	10'	⊙
D	2	8'	⊙
E	1	8'	⊙
F	1	8'	⊙

Inner knee

Remember not to apply SILUET W&M on the knee joints.

Apply gel over the area to be worked using SILUET W&M.


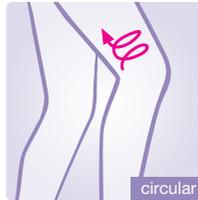
Apply a light circular massage to this area, as shown in the diagram.

Level	Intensity	Time	Mode
A	3	10'	⊙
B	3	10'	⊙
C	2	8'	⊙
D	2	8'	⊙
E	1	8'	⊙
F	1	8'	⊙

Arms

Remember, never apply SILUET W&M on the hands, forearms or joints.

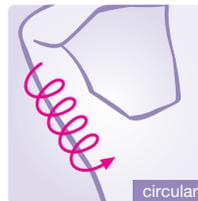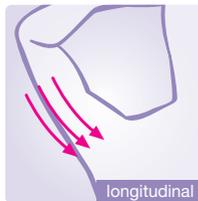Make sure to put a generous amount of gel on the area to be worked, this will improve the effectiveness of the treatment and prevent damage to the machine.

Perform the massage, sliding the handset slowly and without pressing. Massaging this area should be applied toward the body, as shown in the diagram.

Level	Intensity	Time	Mode
A	3	8'	⊙
B	3	8'	⊙
C	2	5'	⊙
D	2	3'	⊙
E	1	3'	⊙
F	1	5'	⊙

CHAPTER 4

Correct, safe use of Siluet W&M

For correct and safe use of SILUET W&M, you should keep the following application times in mind:

- **PER WORK AREA: 15 minutes maximum.**

By “work area” we mean each one of the body parts where SILUET W&M can be applied. The various work areas are: arms, abdomen, waist, hips, inner thighs and inner knees.

E.g.: In a single session you can apply SILUET W&M for 15 minutes on your sides and 15 minutes on your arms. But, for example, you must NEVER perform a 30 minute application session just on your sides in the same day.

- **PER SESSION: 30 minutes maximum.**

By session we mean using SILUET W&M on 2 different work areas (15 minutes per area) in the same day.

- **BETWEEN SESSIONS: 72 hours minimum.**

By “between sessions” we mean the rest time between each session. *E.g.: In one week the session schedule (maximum) would be: Monday, Thursday and Sunday.*

- **EVERY 20 SESSIONS: You should rest for 3 weeks.** After this period you may decide to repeat the first treatment or start the maintenance program.
- **MAINTENANCE SESSIONS:** This depends on the individual, but one session per month should suffice.

VERY IMPORTANT

***** Never exceed the application times (15 minutes per area, 30 minutes total in one day).** Your body does not know “where” you are eliminating the fat from, so a session must never include 15 minutes on your sides, 15 minutes on your legs, 15 minutes on your arms, etc. **REMEMBER:** the maximum time during one day, adding all the parts of the body being treated, must be no more than 30 minutes.

***** And, above all, respect the rest period of 72 hours between sessions.** This is the time that your body needs to eliminate all the fat dissolved during the session through the lymphatic system.

CHAPTER 5

Frequently asked questions

What is cavitation?

Cavitation is a non-invasive (surgery free) body shaping technique that reduces volume by eliminating cellulite and localised fat deposits in the abdomen, hips, flanks, arms, inner thighs and knees.

What are the advantages over other methods?

The main advantage is that we can choose the specific areas of the body that we want to work. It could be said that it is a technology which allows us to concentrate on cellulite selectively.

Is cavitation painful?

Use of SILUET W&M should never cause pain.

Is it normal to hear a slight buzzing in your ears while using SILUET W&M?

Yes, some users may detect a slight buzzing sound during its application.

Can it be used on all parts of the body?

It cannot be applied to all parts of the body. Under no circumstances can it be applied to: the head or face, neck, breasts, the thorax, navel, forearm, joints, genitals, back, ankles or lymph glands.

Does SILUET W&M reduce weight and volume?

SILUET W&M eliminates clusters of fat, which have weight, consequently your body loses volume.

Is SILUET W&M effective on orange peel syndrome?

Yes, SILUET W&M is effective on orange peel syndrome because it destroys adipose cells that cause fat storage.

Can SILUET W&M be used to remove wrinkles?

No, SILUET W&M has not been designed to remove wrinkles.

Can SILUET W&M be used on the face?

SILUET W&M must NEVER be used on the face, eyes or neck.

After how many sessions will I begin to see results?

You should notice results after the first five sessions, however results may vary depending on the subject.

How long should I wait between sessions?

It is advisable to wait 72 hours between sessions.

Does cellulite reappear after using SILUET W&M?

Although SILUET W&M ensures effective and long-lasting results, we cannot guarantee against cellulite reappearing. To prevent the appearance of new cellulite we recommend a balanced diet, regular physical exercise and maintenance sessions using SILUET W&M.

Are there any side effects?

If you are healthy, SILUET W&M has no reported side effects, nevertheless in some people it may produce a reddening of the skin that can last for about 2 hours.

Is it necessary to combine SILUET W&M with a special diet?

No, it is not necessary to follow any special diet, but a balanced diet and proper hydration is highly advisable.

How much water should I drink?

It is essential to begin hydrating yourself two hours before applying SILUET W&M. The recommended quantity is at least 2 litres during these periods, as it helps to eliminate toxins.

Is it necessary to combine SILUET W&M with a sports activity?

No, this is not strictly necessary, but at least 30 minutes of physical activity will enhance the results as it helps to drain fluids.

Which sports activities can supplement SILUET W&M?

The vibration platform is the ideal complement due to its toning and draining qualities. But if you do not have access to a platform, then walking, exercise bike, swimming or any other aerobic exercise for at least 30 minutes is very helpful.

How can I enhance the effects of SILUET W&M?

By combining a balanced diet with regular physical exercise and drinking plenty of water. Added to which, drainage massage can also be very useful.

What are the anti-cellulite and toning gels for? When do I use them?

1. Anti-cellulite gel: Contains active ingredients that help eliminate cellulite.
2. Toning gel: Contains active ingredients that make the treated area firmer.

The application sequence is as follows: during the first few days of treatment use the anti-cellulite gel and once you begin to notice results, use the toning gel. This is applied before the Cavitation gel, as with any other serum. The area should be massaged until it is well absorbed and then apply plenty of cavitation gel to start the session.

Can I use SILUET W&M without applying gel?

NEVER, because the gel acts as a conductor for the ultrasonic sound waves with your body. If you do not use gel then the SILUET W&M application loses its effectiveness and the machine can be severely damaged.

Can I use any type of gel?

We recommend using only COSMOBELLE gel, as the use of any other type may damage the machine.

When is it inadvisable to use SILUET W&M?

- If you are pregnant or think you might be.
- If you are breast-feeding.
- Undergoing fertility treatments.
- If you have any kind of heart disease.
- If you suffer from skin problems.
- Osteoporosis.
- Fibromyalgia.
- Malignant tumours.
- Hypertension.
- High cholesterol levels.
- Kidney problems.
- Physical shortcomings.
- Metal implants or prostheses.
- Hearing implants.
- Physical malformations.
- Epilepsy.
- Long-term treatment with medications: anti-inflammatories, anticoagulants, antibiotics or antihistamines.
- Depression.
- Psychiatric or psychological treatment.
- Liver problems (fatty liver).
- Irritable bowel syndrome.
- Diabetes.

Can I use SILUET W&M if I am pregnant?

NEVER use SILUET W&M if you are pregnant or suspect you might be.

Can SILUET W&M be used during lactation?

Do not use it during lactation.

Can SILUET W&M be used on the breasts?

NEVER apply it to the breasts.

Can I use SILUET W&M during menstruation?

Use on the abdominal area is not advisable during menstruation as it may cause increased bleeding. However, it can be applied in other areas.

Can I use SILUET W&M if I have an IUD?

Yes.

Can SILUET W&M be used on tattooed skin?

No. You should avoid passing the handset over tattooed areas.

Can I use SILUET W&M if I have implants (pins, etc)?

SILUET W&M can be used if you have implants as long as you NEVER massage these areas: for example, if the implant is in the knee, there is no problem massaging the abdomen. Nor will there be problem if the implant is in your back, because you should never use SILUET W&M on this part of the body.

Can SILUET W&M affect my cholesterol levels?

Use of SILUET W&M may increase your cholesterol levels.

Can SILUET W&M be used with skin conditions?

It should not be used if you have any skin condition. Consult your doctor.

Can I use SILUET W&M if I have hearing implants?

Do not use it if you have hearing implants fitted.

Can I use SILUET W&M when I am bathing or showering?

Under no circumstances should the handset or the machine come into contact with water due to the risk of electrocution.

How do I clean SILUET W&M and the handset?

Never immerse SILUET W&M or the handset in water to clean it. Use a cloth slightly dampened with water.

CHAPTER 6

Important safety advice

Do not modify the electric plug supplied with the unit. If it does not fit the wall socket then get an electrician to install a suitable electrical socket.

This product should be used on a 220-240 volt circuit and has an earthed plug similar to that shown in **fig. 5.1**

Follow these recommendations:

1. Never use accessories that are not recommended by the manufacturer.
2. Do not use the machine if the mains cable becomes damaged or worn.
3. Keep the mains cable away from hot surfaces.
4. Do not use the machine outdoors.
5. To unplug the machine, first put the switch to the «0» position and then unplug it from the wall socket.




Protection

This model is equipped with an electronic circuit breaker to protect both the electronic and electrical circuits.

This safety device (fuse) will be automatically triggered if the unit is subjected to abnormal force or heat.

Should this occur:

1. Unplug the cable from the mains supply.
2. Set switch (L), located on the side of the machine, to the «0» (OFF) position **fig. 5.2**
3. Unscrew the machine's outer casing
4. Change the fuse (1 Amp) located inside the fuse holder, next to the mains cord connector **fig. 5.3**
5. Plug the unit into a wall socket. Set switch (L) to the «I» (ON) position **fig. 5.2**

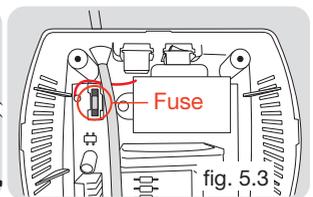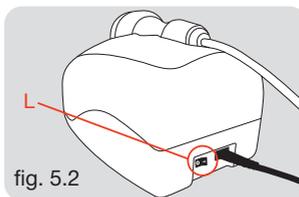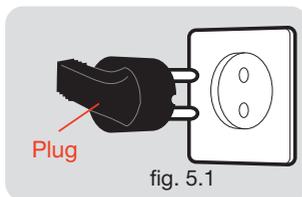Warranty

We suggest that you read through this document carefully in order to fully understand the cover that the manufacturer offers as well as your responsibilities as owner.

Guarantees

GENERAL INFORMATION

EXERCYCLE, S.L., as distributor of Copsmobelle, would like to thank you for purchasing this product and the trust you have placed in us. Our products have been subjected to rigorous quality checks. Nevertheless, should you experience any problem with this appliance during the guarantee period then please read through these warranty conditions and then get in touch with your local authorised dealer.

This guarantee certificate has been drafted in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25th May, 1999, regarding on certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10th July; Guarantees on the sale of consumer goods.

To this effect, EXERCYCLE, S.L. VAT n° B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product under the following terms and conditions:

1. GUARANTEE PERIOD

The product covered by this guarantee certificate is protected for a period of TWO (2) YEARS as from the date of receiving the product, which will be taken as the date that appears on the sales receipt for the product or on the corresponding delivery docket should this be later. It is essential that this be presented in the event of possible repairs. The policy regarding the warranty period may vary between different countries. The stated warranty period is only applicable to Europe.

2. RIGHTS UNDER THE GUARANTEE

- Repair or replacement of the product which, in such instance, will be carried out within a reasonable period and without significant inconvenience to the consumer, taking into account of the nature of the goods and the purpose for which the consumer required them.
- Those cases in which the repair or replacement represents a disproportionate cost for the seller in comparison with other forms of remedy will be exempt assuming these do not prejudice the consumer.
- Repairs will be free for the consumer with regard to labour, materials and carriage.
- In the event that the machine is repaired or replaced with a new one, the guarantee period will be suspended temporarily until the new appliance has been delivered and will carry at least six months of guarantee, however under no circumstances will a new 2 year period be initiated.

- The buyer may, as is apparent, require (i) a reduction in price or proceed to a (ii) recession of the contract when the repair or replacement is not carried out within a reasonable time or with significant inconvenience to the buyer. The reduction in price (i) will be proportional to the difference between the value that the goods would have had if they had been delivered in accordance with the contract and the actual value that the goods were worth at the time of the said delivery. Recession of the contract (ii) will not proceed when the non-conformity is of minor importance.

3. THE GUARANTEE DOES NOT COVER

- Defects attributable to use other than that established in the User Manual.
- Professional use or commercial leasing or hiring of the product.
- Poor installation, poor assembly or breakage of components due to improper installation, (the User Manual gives detailed instructions on how to install and use the machine).
- Improper use of the product (use without applying gel, exceeding operating times, exposure to high moisture levels, etc.) or problems arising from poor maintenance (accumulation of dirt on any part of the product or its accessories, etc.) (check the User Manual for detailed instructions on how to install and use the product).
- Repairs carried out by persons not authorised by EXERCYCLE, S.L. or appliances that have been opened or tampered with by persons not authorised by the manufacturer.
- Deviations from the standard product, such as light vibration, noise or breakage of components that do not affect the performance of the machine and which cannot be attributed to manufacturing faults.
- Scratches, scrapes, abrasions, knocks, falls, etc. caused by carelessness.
- Defects on the appliance caused by damage during transit.

4. CLAIMING AGAINST THE GUARANTEE

- Claims made for nonconformities will be limited to three years as from the date that the goods were delivered.
- The consumer must report the existence of the defect within a period of 2 months as from the date that the defect was discovered.
- This guarantee does not affect the consumer's statutory rights under any applicable law, unless it is incompatible with the application of other actions derived from compensation sought for hidden faults in the contract of sale.
- With regard to queries that may arise from the interpretation or application of this guarantee, these will be governed by the corresponding national law.

Any discrepancy arising from the interpretation, validity or application of this contract will be subject to the courts and tribunals of Vitoria, Spain.

Certificates

By virtue of this document, Exercycle S.L., as distributor for Cosmobelle, declares that this product meets the essential requirements and other relevant provisions of Directives 2004/108/EC and 2006/95/EC.


Parts list

En

2 Handset

16 Mains cable

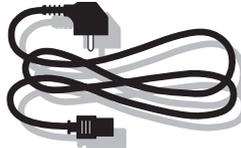G01 Main unit

To order replacement parts, please indicate:

The machine model	Corresponding part N°	Quantity
-------------------	-----------------------	----------

Ejemplo:

CB1500	16	1
--------	----	---

Technical specifications

Net weight: 2.3 kg.

Dimensions (length x width x height): 32 x 17 x 13 cm.

Freq. Ultrasound: 40 kHz

Output power: 5 to 20 W (depending on program)

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