Cavislim

cavitation

Tecnovita
by **BH**

ENGLISH

FRANÇAIS

cavitation safe and effective figure shaping
modelage sûr et efficace de la silhouette

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Congratulations. First we would like to thank you for purchasing CAVISLIM and for the trust you have placed in us. By means of this handy User Manual and following its simple instructions you will achieve results comparable to those obtained from the hands of a professional beautician. We strongly recommend that you read this User Manual thoroughly. Reading and understanding it will help you use CAVISLIM effectively and safely.

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content

chapter 1

- contents of the box
- assembling Cavislim

chapter 2

- basic characteristics of Cavislim
- cavislim operating instructions

chapter 3

- cavislim application modes

chapter 4

- for correct, safe use of Cavislim

chapter 5

- frequently asked questions

chapter 6


- important safety instructions
- warranty
- certificates
- technical specifications
- parts list

chapter 1


contents of the box

Inside this box you will find everything you need to enjoy CAVISLIM from now onward. Carefully check the following illustration showing all of the elements that the pack should contain.

Hand unit support.

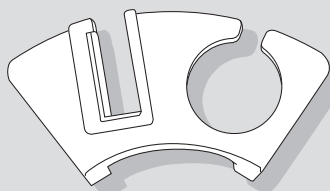Hand unit

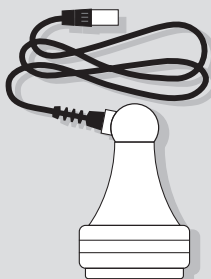
User manual


DVD


Main unit

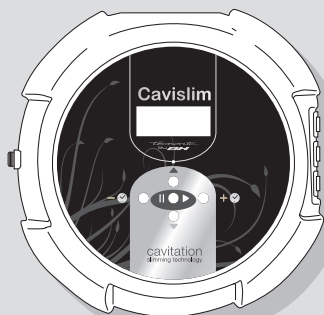Operating pedal

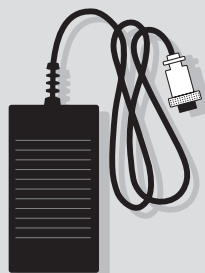Cavitation gel

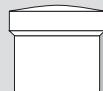
Toning serum


Allen key

4 screws

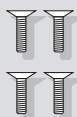Mains cable

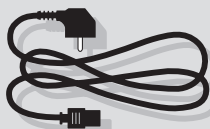
Anti-cellulite serum

assembling Cavislim

First we are going to fit the **hand unit support** onto the **main unit**. To do this, place the main unit on a table and position the **hand unit support** by inserting the **screws** as shown in the diagram.

Using the **Allen key** tighten the **screws** securely so that the two parts fit together perfectly.

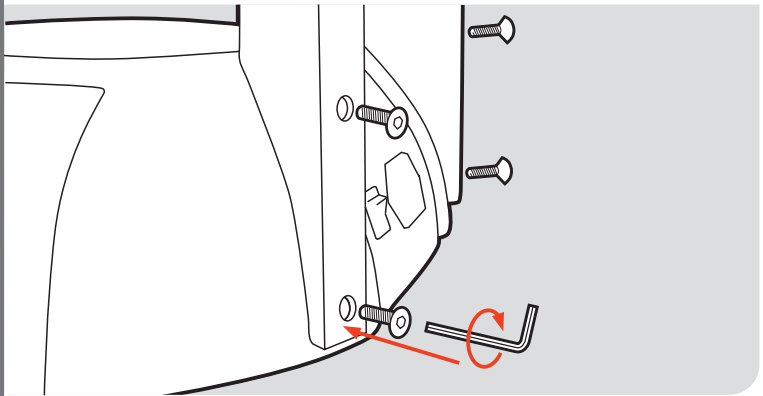basic characteristics of Cavislim

CAVISLIM and CAVISLIM EXCEL are a high-tech devices that allow us to selectively and efficiently remove the build-up of fat (cellulite, orange skin syndrome) from our body. The technology used is that of ultrasound emission which produces the physical phenomenon known as cavitation, causing the removal of fat cells. The technology with proven results is widely used in beauty salons all over the world in order to fight cellulite painlessly and non-invasively.

Characteristics for model CAVISLIM Ref. YS30

Time: Programmable from 1 to 15 minutes.

Intensity: 4 intensity levels (5-10-15-20)

Mode: 4 programmable pulsed work modes (P1-P2-P3-P4)

Characteristics for model CAVISLIM EXCEL Ref. YS50

Time: Programmable from 1 to 15 minutes.

Intensity: 5 intensity levels (5-10-15-20-25)

Mode: 4 programmable pulsed work modes (P1-P2-P3-P4)

1 continuous work mode (P5)

With these units you only have to worry about programming in three simple functions: intensity, time and work mode. The intensity function controls how deep the sound waves penetrate your skin. The greater the intensity, the deeper the waves reach. The time function allows you to choose how long you wish to work with CAVISLIM. For reasons of safety the time function allows you to choose up to 15 minutes which corresponds to the maximum recommended time for working a specific area. However, you can reduce this maximum time by 1 minute intervals. Once the selected time has elapsed CAVISLIM will alert you with a beep and a message on the display showing the text "Time out".




CAVISLIM also lets you choose the work mode, helping you select the time interval during which the waves reach your skin.

Cavislim operating instructions

Precautions: Consult your doctor before beginning to use this appliance. This warning is especially important for people with health problems. Use by those under 18 years of age is not recommended. Keep this appliance out of children's reach.

Start up

- Connect the hand unit to the main unit. To do this, take the plug for the **hand unit** and insert into the socket on the main body as shown in Fig. 2.1.
- Next, take the jack for the **operating pedal** and insert into the socket on the **main body** as shown in Fig. 2.2. Check that all the connections are made correctly before connecting it to the mains supply.
- Plug one end of the mains cable into the socket on the machine and the other end into the electric wall socket (220V) Fig. 2.3. Switch the machine to position (I) to start the machine.

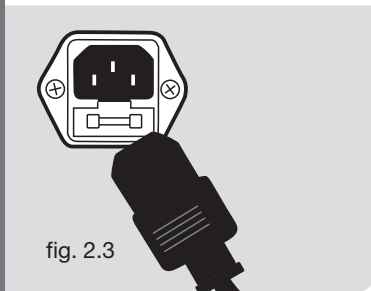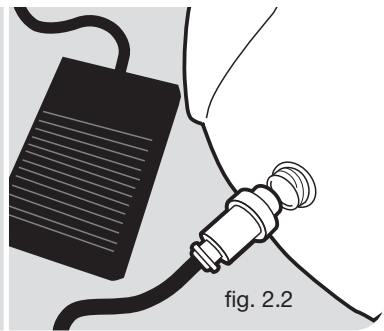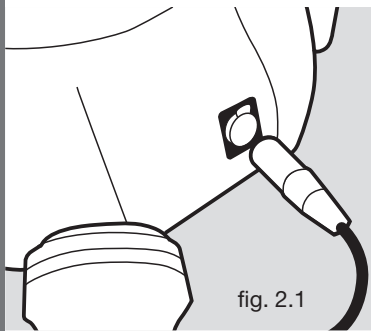
Monitor elements

The CAVISLIM monitor consists of a display screen, an LED and five buttons or keys.

- The display provides information about the cavitation session.
- The ON LED (green) tells you that the device is operating.
- The keys have the following functions:


- ▲ and ▼ change the work **MODE**
- + and - adjust **TIME** and **INTENSITY**
- || pause and ▶ start session

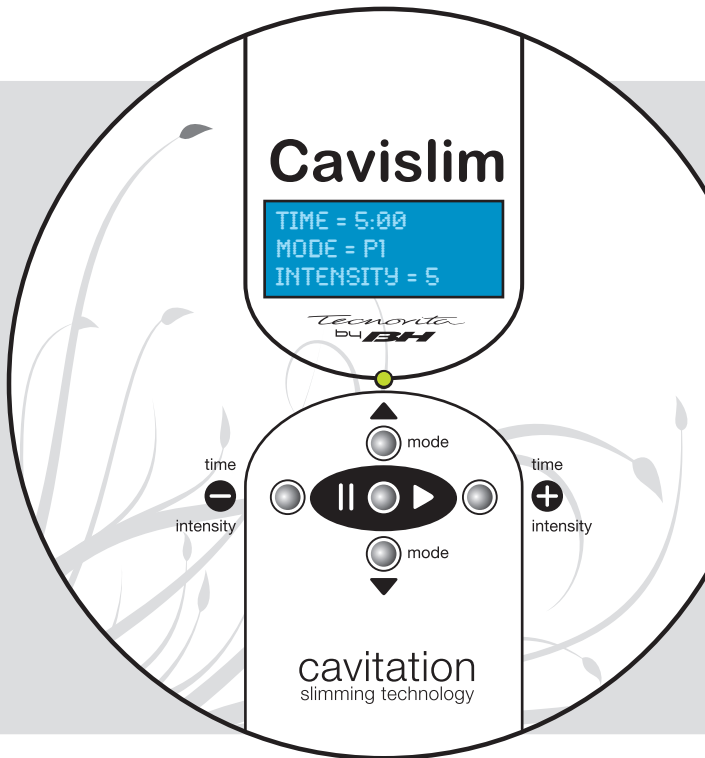Cavitation session

There are 3 variables for the massage:

- **Time:** measures the session in minutes and seconds. Choose between 1 to 15 minutes (in steps of 1 minute).
- **Intensity:** establishes the intensity level for the session. On Mod. **YS30** you can select from levels 5, 10, 15 and 20 whereas with Mod. **YS50** you can select from levels 5, 10, 15, 20 and 25.
- **Work mode:** Pulse type (impulse emission frequency). Mod. **YS30** allows you to select from levels P1, P2, P3 and P4. Mod. **YS50** allows you to select from levels P1, P2, P3, P4 and P5. Impulse emission frequency: P1 very low
P2 low
P3 high
P4 very high
P5 continuous

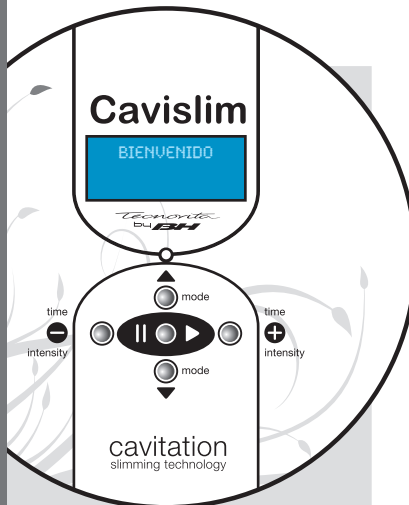Beginning a session

Once you have plugged the machine into the mains supply and switched it to position (I) the following message will appear on the monitor: **BIENVENIDO (WELCOME)**

This message will remain visible for 3 seconds, after which the following default settings will appear:

TIME: 10:00
MODE: P1
INTENSITY: 5

From here you can modify the time and work mode settings to establish the type of session you want:

Changing session duration

Once you are on the default setting screen you can modify the duration of the session using



+ and **-**


The displayed default value is 10 minutes but this can be changed in intervals of 1 minute between the values 1 and 15.

Changing work mode

On the same default screen you can change the work mode for your session using the **▲** and **▼** keys.

The default value shown is P1 (very low pulse frequency) but this can be modified by selecting from the values **P1**, **P2**, **P3** and **P4** (model YS30) and **P1**, **P2**, **P3**, **P4** and **P5** (only on model YS50).

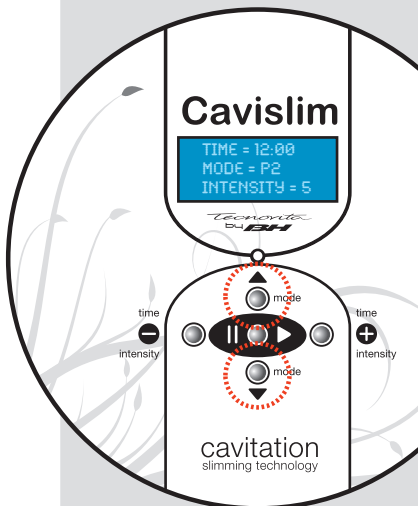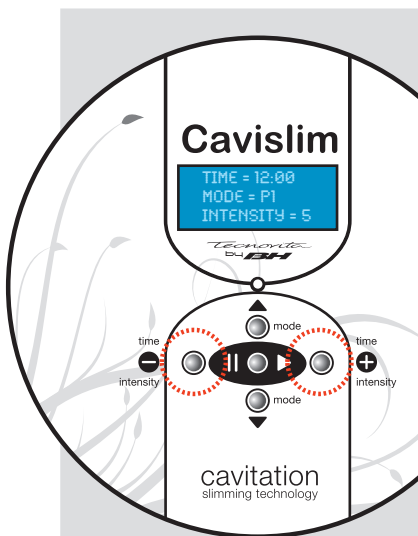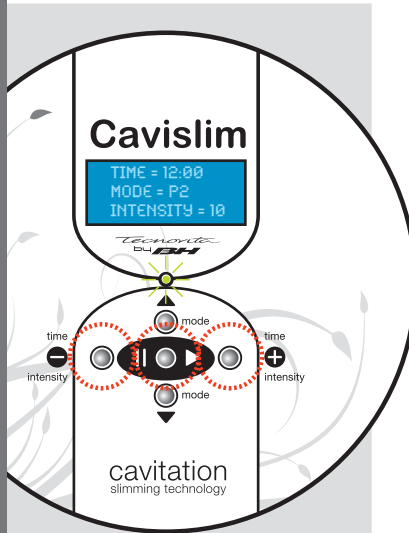Changing message intensity

Once you have selected the desired time and work mode settings, you can set the last of the parameters: Intensity.

To do this, press the button  and with the machine on (a green light appears in the centre of the unit), you can use the **+** and **-** keys to select the intensity level.

The default intensity level is 5 but this can be changed by selecting from 5, 10, 15 and 20 (model YS30) and 5, 10, 15, 20 and 25 (only on model YS50).

After pressing the key  it is not possible to modify the time setting.

Finally, take hold of the hand unit and push down on the pedal to begin your cavitation session.

When you tread on the pedal you will see the time countdown on the screen and the hand unit will begin to emit impulses. Place the hand unit on the area you wish to work and massage it as shown on the DVD in the “Application modes” section.


If you wish to pause the session, just stop pressing on the pedal and the impulse emission will stop immediately.

When the time counter reaches 0.00 the unit will stop and emit a series of beeps to tell you that the selected session has finished. The following message will appear on the monitor:

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Pressing any key on the machine will take you back to the default screen.

TIME: 10:00
MODE: P1
INTENSITY: 5

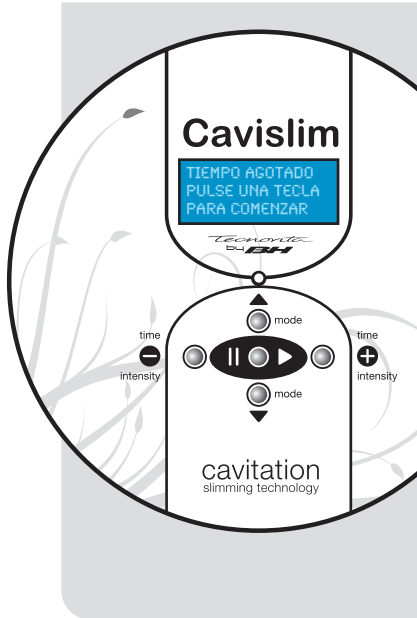chapter 3

Cavislim application modes

Always apply an abundant amount of gel to the area and begin to work it by sliding the hand unit slowly and without pressing, as shown in the diagrams.


If you need to purchase gel or any of our preparation creams you can find them in the shop where you purchased CAVISLIM.

In each of the following sections you will find a table with application protocols. You should decide which level applies to you.

- A** Large accumulation of high density fat
- B** Large accumulation of medium density fat
- C** Medium accumulation of medium density fat
- D** Small accumulation of fat
- E** orange peel syndrome
- F** Maintenance sessions

abdomen central


Remember, do not to apply CAVISLIM to the navel, breast or the thorax area

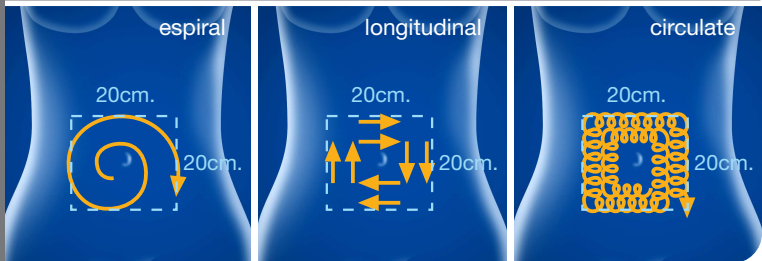Slide the hand unit slowly, without pressing. For best results we recommend applying CAVISLIM by combining different massage techniques (spiral, crossways, circular) in different sessions.

Nivel	YS50		YS30		Time
	Mode	Intensity	Mode	Intensity	
A	3	10	3	10	10'
B	3	10	3	10	10'
C	3	10	3	10	7'
D	2	10	2	10	7'
E	2	10	2	10	7'
F	2	10	2	10	7'

sides

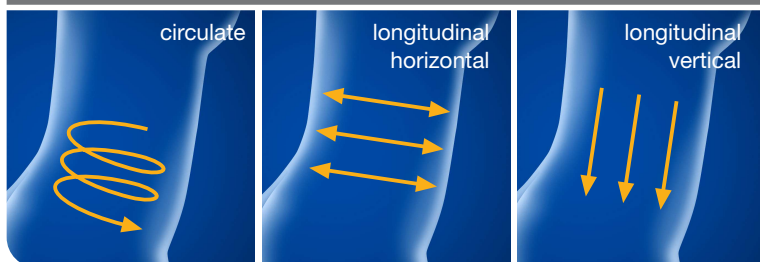
Circular and longitudinal massage should be applied from top to bottom as shown in the diagrams. It is important to massage in the correct direction as this will promote drainage.

Nivel	YS50		YS30		Time
	Mode	Intensity	Mode	Intensity	
A	4	15	4	15	10'
B	3	10	3	10	10'
C	2	10	2	10	7'
D	1	5	1	5	7'
E	1	5	1	5	5'
F	2	5	2	5	5'

Apply gel over the area to be worked using CAVISLIM

hips

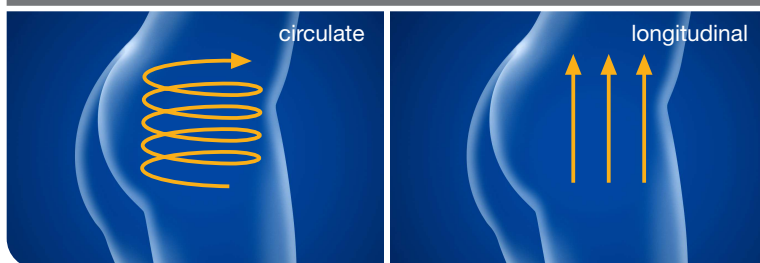
Both the circular and longitudinal massage should be applied from bottom to top as shown in diagrams. This will produce an effect on a larger accumulation of fat during treatment.

Nivel	YS50		YS30		Time
	Mode	Intensity	Mode	Intensity	
A	3	10	3	10	7'
B	2	10	2	10	7'
C	2	5	2	5	7'
D	1	5	1	5	5'
E	1	5	1	5	5'
F	1	5	1	5	5'

Apply gel over the area to be worked using CAVISLIM

buttocks

Apply gel over the area to be worked using CAVISLIM

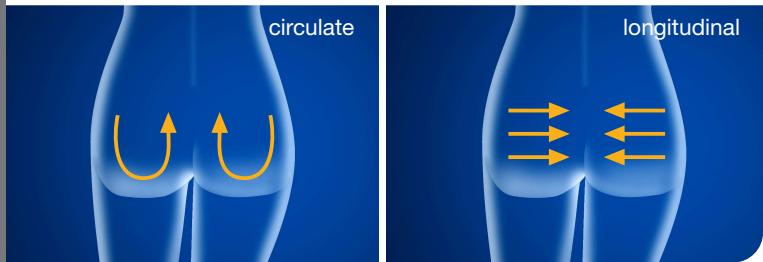
Longitudinal and circular massage should be done in the direction of the arrows, as shown in the diagrams.

Nivel	YS50		YS30		Time
	Mode	Intensity	Mode	Intensity	
A	5	25	4	20	10'
B	5	15	4	15	7'
C	4	20	4	20	7'
D	3	15	3	15	7'
E	2	15	2	15	10'
F	3	15	3	15	10'

inner thigh

Remember, do not to apply CAVISLIM from the knee to the foot, there is no significant accumulation of fat in this area.

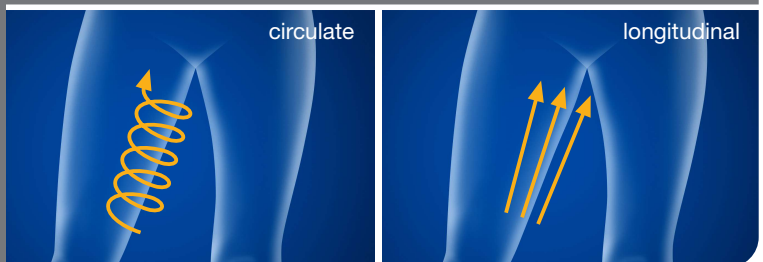
Both the longitudinal and circular massage should be done from bottom to top, as shown in the diagrams. Follow these simple instructions to help you get the best results.

Nivel	YS50		YS30		Time
	Mode	Intensity	Mode	Intensity	
A	3	10	3	10	10'
B	3	10	3	10	10'
C	3	10	3	10	7'
D	2	10	2	10	7'
E	2	10	2	10	7'
F	2	10	2	10	7'

back of the knee

EN

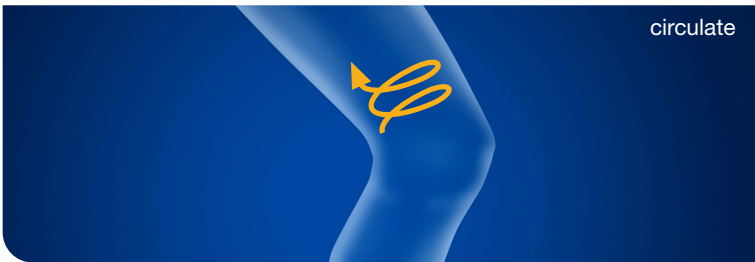
Apply a light circular massage to this area, as shown in the diagram.

Nivel	YS50		YS30		Time
	Mode	Intensity	Mode	Intensity	
A	3	5	3	5	10'
B	3	5	3	5	10'
C	3	5	3	5	7'
D	2	5	2	5	7'
E	2	5	2	5	7'
F	2	5	2	5	7'

Remember, do not apply CAVISLIM on the kneecap. Apply gel on the area to be worked using CAVISLIM

arms

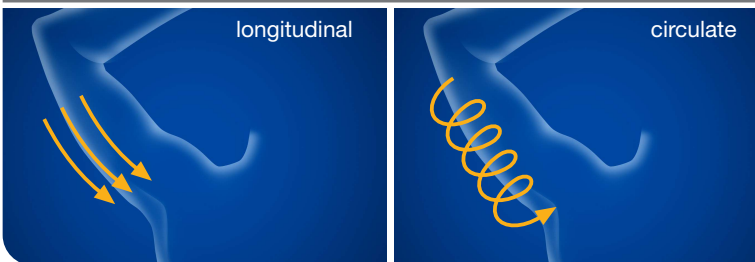Put a generous amount of gel on the area to be worked and this will improve the effectiveness of the treatment. Massaging this area should be applied toward the body, as shown in the diagram.

Nivel	YS50		YS30		Time
	Mode	Intensity	Mode	Intensity	
A	2	5	2	5	7'
B	2	5	2	5	5'
C	1	5	1	5	5'
D	1	5	1	5	3'
E	1	5	1	5	3'
F	1	5	1	5	5'

Remember, never apply CAVISLIM on the hands, forearms or joints.

chapter 4

for correct and safe use of Cavislim

For correct and safe use of Cavislim, you should keep the following application times in mind:

- **Per work area:**

15 minutes maximum.

By "work area" we mean each one of the body parts where CAVISLIM can be applied. The various work areas are: arms, abdomen, waist, hips, inner thighs and inner knees.

E.g.: In a single session you can apply Cavislim for 15 minutes on your sides and 15 minutes on your arms. But, for example, you must NEVER perform a 30 minute application session just on your sides in the same day.

- **Per session:**

30 minutes maximum.

By session we mean using CAVISLIM on 2 different work areas (15 minutes per area) in the same day.

- **Between sessions:**

72 hours minimum.

By "between sessions" we mean the rest time between each session.

E.g.: In one week the session schedule (maximum) would be: Monday, Thursday and Sunday.

- **Every 20 sessions:**

You should rest for 3 weeks.

After this period you may decide to repeat the first treatment or start the maintenance program.

- **Maintenance sessions:**

This depends on the individual, but one session per month should suffice.

chapter 5

frequently asked questions

What is cavitation?

Cavitation is a non-invasive (surgery free) body shaping technique that reduces volume by eliminating cellulite and localised fat deposits in the abdomen, hips, flanks, arms, inner thighs and knees.

What are the advantages over other methods?

The main advantage is that we can choose the specific areas of the body that we want to work. It could be said that it is a technology which allows us to concentrate on cellulite selectively.

Is cavitation painful?

Use of CAVISLIM should never cause pain.

Is it normal to hear a slight buzzing in your ears while using CAVISLIM?

Yes, some users may detect a slight buzzing sound during its application.

Can it be used on all parts of the body?

It cannot be applied to all parts of the body. Under no circumstances can it be applied to the head or face, neck, breasts, the thorax, navel, forearm, joints, genitals, back, ankles or lymph glands.

Does CAVISLIM reduce weight and volume?

CAVISLIM eliminates clusters of fat, which have weight, consequently your body loses volume.

Is CAVISLIM effective on orange peel syndrome?

Yes, CAVISLIM is effective on orange peel syndrome because it destroys adipose cells that cause fat storage.

Can CAVISLIM be used to remove wrinkles?

No, CAVISLIM has not been designed to remove wrinkles.

Can CAVISLIM be used on the face?

CAVISLIM must never be used on the face or eyes.

After how many sessions will I begin to see results?

You should notice results after the first five sessions, however results may vary depending on the individual.

How long should I wait between sessions?

It is advisable to wait 72 hours between sessions.

Does cellulite reappear after using CAVISLIM?

Although CAVISLIM ensures effective and long-lasting results, we can not guarantee against cellulite reappearing. To prevent the reappearance of

cellulite we recommend a balanced diet, regular physical exercise and maintenance sessions using CAVISLIM.

Are there any side effects?

If you are healthy, CAVISLIM has no reported side effects, nevertheless in some people it may produce a reddening of the skin that can last for about 2 hours.

What are the anti-cellulite and toning gels for? When do I use them?

1. Anti-cellulite gel: Contains active ingredients that help eliminate cellulite.
2. Toning gel: Contains active ingredients that make the treated area firmer.

The application sequence is as follows: during the first days of application we use the anti-cellulite gel and the toning gel during the last days of application. This is applied before the Cavitation gel. It is applied like any other serum. Make sure it is fully absorbed. Then apply the cavitation gel to begin the session.

Is it necessary to combine CAVISLIM with a special diet?

No, it is not necessary to follow any special diet, but a balanced diet and proper hydration is highly advisable.

How much water should I drink?

It is essential to begin hydrating yourself two hours before applying CAVISLIM. It is also advisable to drink more water after finishing a session. The recommended quantity is at least 2 litres during these periods, as it helps to eliminate toxins.

Is it necessary to combine CAVISLIM with a sports activity?

No, this is not strictly necessary, but at least 30 minutes of physical activity will enhance the results as it helps to drain fluids.

Which sports activities can supplement CAVISLIM?

The vibration platform is the ideal complement due to its toning and draining qualities. But if you do not have access to a platform, then walking, exercise bike, swimming or any other aerobic exercise for at least 30 minutes is very helpful.

How can I enhance the effects of CAVISLIM?

By combining a balanced diet with regular physical exercise. Added to which, drainage massage can also be very useful.

Can I use CAVISLIM without applying gel?

Never, because the gel acts as a conductor for the ultrasonic sound waves with your body. If you do not use gel then the CAVISLIM application loses its effectiveness and the machine can be severely damaged.

Can I use any type of gel?

We recommend using only TECNOVITA by BH gel, as the use of any other type may damage the machine.

When is it inadvisable to use CAVISLIM?

- if you are pregnant or think you might be
- during lactation
- undergoing fertility treatments
- if you have any kind of heart disease
- if you suffer from skin problems
- osteoporosis
- fibromyalgia
- malignant tumours
- hypertension
- high cholesterol levels
- diabetes
- kidney problems
- physical shortcomings
- metal implants or prostheses.
- physical malformations.
- epilepsy
- long-term treatment with medications: anti-inflammatories, anticoagulants, antibiotics or antihistamines
- depression
- psychiatric or psychological treatment
- liver problems (fatty liver)
- irritable bowel syndrome

Can I use CAVISLIM if I am pregnant?

No, you must never use CAVISLIM if you are pregnant or think you might be.

Can CAVISLIM be used during lactation?

Do not use it during lactation.

Can CAVISLIM be used on the breasts?

It must never be applied to the breasts.

Can I use CAVISLIM during menstruation?

Use on the abdominal area is not advisable during menstruation as it may cause increased bleeding. However, it can be applied in other areas.

Can I use CAVISLIM if I have an IUD?

Yes

Can CAVISLIM affect my cholesterol levels?

Use of CAVISLIM may increase your cholesterol levels.

Can CAVISLIM be used with skin conditions?

It should not be used if you have any skin condition. Consult your doctor.

Can I use CAVISLIM if I have hearing implants?

Do not use it if you have implants fitted.

Can I use CAVISLIM when I am bathing or showering?

Under no circumstances should the handset or the machine come into contact with water due to the risk of electrocution.

How do I clean CAVISLIM and the handset?

Never immerse CAVISLIM or the handset in water. To clean it, use a cloth slightly dampened with water.

chapter 6

important safety instruction

This product is supplied with a mains cable that has an earth wire and a plug with an earth pin. In the event that the appliance malfunctions then the earth connection will redirect the electrical current away from the machine and reduce the risk of an electrical shock.

DANGER! If the equipment is not earthed correctly then there is a risk of electrical shock.

Do not modify the electric plug supplied with the unit. If it does not fit the wall socket then get an electrician to install an electrical socket that does.

This product should be used on a 220-240 volt circuit and has an earthed plug similar to that shown in Fig 5.1 Follow these instructions:

1. Never use accessories that are not recommended by the manufacturer.
2. Do not use the machine if the mains cable becomes damaged or worn.
3. Keep the mains cable away from hot surfaces.
4. Do not use the machine outdoors.
5. To unplug the machine, first put the switch to the «0» position and then unplug it from the wall socket.

Protection

This model is equipped with an electronic circuit breaker to protect both the electronic and electrical circuits.

This safety device (fuse) will be automatically triggered if the unit is subjected to abnormal force or heat.

Should this occur:

1. Unplug the cable from the mains supply.
2. Set switch (L), located on the casing, to the «0» (OFF) position. Fig. 5.2.
3. Change the fuse (C) (1 Amp).
4. Plug the unit into a wall socket equipped with an EARTH terminal. Set switch (L) to the «I» Position, Fig. 5.2

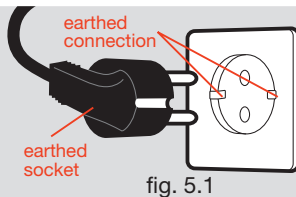
fig. 5.1

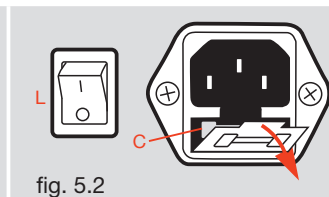
fig. 5.2

guarantee

We suggest that you read through this document carefully in order to fully understand the cover that the manufacturer offers with this guarantee as well as your responsibilities as owner.

GUARANTEES

GENERAL INFORMATION

EXERCYCLE, S.L. would like to thank you for purchasing this product and the trust you have placed in us. Our products have been subjected to rigorous quality checks. Nevertheless, should you experience any problem with this appliance during the guarantee period then please read through these warranty conditions and then get in touch with your local authorised dealer. This guarantee certificate has been drafted in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25th May, 1999, regarding on certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10th July; Guarantees on the sale of consumer goods. To this effect, EXERCYCLE, S.L. VAT nºbc B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product under the following terms and conditions:

1. GUARANTEE PERIOD

The product covered by this guarantee certificate is protected for a period of TWO (2) YEARS as from the date of receiving the product, which will be taken as the date that appears on the sales receipt for the product or on the corresponding delivery docket should this be later. It is essential that this be presented in the event of possible repairs. The policy regarding the warranty period may vary between different countries. The stated warranty period is only applicable to Europe.

2. RIGHTS UNDER THE GUARANTEE

- Repair or replacement of the product which, in such instance, will be carried out within a reasonable period and without significant inconvenience to the consumer, taking into account of the nature of the goods and the purpose for which the consumer required them.
- Those cases in which the repair or replacement represents a disproportionate cost for the seller in comparison with other forms of remedy will be exempt assuming these do not prejudice the consumer.
- Repairs will be free for the consumer with regard to labour materials and postage.
- In the event that the machine is repaired or replaced with a new one, the guarantee period will be suspended temporarily until the new appliance has been delivered and will carry at least six months of guarantee, however under no circumstances will a new 2 year period be initiated.
- The buyer may, as is apparent, require (i) a reduction in price or proceed to a (ii) recession of the contract when the repair or replacement is not carried out within a reasonable time or with significant inconvenience to the buyer.

The reduction in price (i) will be proportional to the difference between the value that the goods would have had if they had been delivered in accordance with the contract and the actual value that the goods were worth at the time of the said delivery. Recession of the contract (ii) will not proceed when the non-conformity is of minor importance.

3. THE GUARANTEE DOES NOT COVER

- Defects attributable to use other than that established in the User Manual.
- Professional use or commercial leasing or hiring of the product.
- Poor installation, poor assembly or breakage of components due to improper installation, (the User Manual gives detailed instructions on how to install and use the machine).
- Improper use of the product (use without applying gel, exceeding operating times, exposure to high moisture levels, etc.) or problems arising from poor maintenance (accumulation of dirt on any part of the product or its accessories, etc.) (check the User Manual for detailed instructions on how to install and use the product).
- Repairs carried out by persons not authorised by EXERCYCLE, S.L. or appliances that have been opened or tampered with by persons not authorised by the manufacturer.
- Deviations from the standard product, such as light vibration, noise or breakage of components that do not affect the performance of the machine and which cannot be attributed to manufacturing faults.
- Scratches, scrapes, abrasions, knocks, falls, etc. caused by carelessness.
- Defects on the appliance caused by damage during transit.

4. CLAIMING AGAINST THE GUARANTEE

- Claims made for nonconformities will be limited to three years as from the date that the goods were delivered.
- The consumer must report the existence of the defect within a period of 2 months as from the date that the defect was discovered.
- This guarantee does not affect the consumer's statutory rights under any applicable law, unless it is incompatible with the application of other actions derived from compensation sought for hidden faults in the contract of sale.
- With regard to queries that may arise from the interpretation or application of this guarantee, these will be governed by the corresponding national law. Any discrepancy arising from the interpretation, validity or application of this contract will be subject to the courts and tribunals of Vitoria, Spain.

certificates

By virtue of this document, Exercycle S.L. declares that this product meets the essential requirements and other relevant provisions of Directives 2004/108/EC and 2006/95/EC

parts list

EN


1 Operating pedal

2 Hand unit


3 Hand unit support.

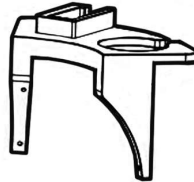16 Mains cable

G01 Main unit

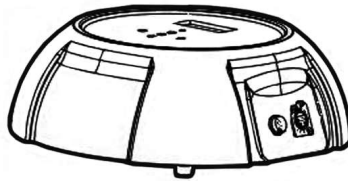To order replacement parts, please indicate:

Model of the machine	Corresponding N° of the piece	Quantity
----------------------	-------------------------------	----------

Example:

YS30	16	1
------	----	---

technical specifications

Weight: 7 kg

Dimensions: L x W x H: 34 x 38 x 13 cm.

Ultrasound: 40 kHz

Output power: 5 to 25 Watts (depending on program)

Power supply: Voltage: AC 220V / Frequency: 50 Hz

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